

FERNO Compact 2 Track



In-vehicle storage bracket

Cabin-mounted two-part bracket holds both the chair and track - so they're to hand when you need them and secure when you don't



**FOR MORE DETAILED
INFORMATION SEE OVER**

Specifications

FERNO – is a world leading innovator and manufacturer of patient transport solutions.

Our product range is born of a closer understanding of our customers' needs and the practical demands of their work. Ferno's industry-leading **Compact 2 Track Bracket** typifies our passion for top quality, highly functional design... patient transfer products that don't just do the job, they do it **BETTER**.

Width	Depth	Height	Total Height Required	Bracket weight
435 mm With chair 450 mm	180 mm With chair 290 mm	687 mm With chair 820mm	850 mm for safe removal	7.5kg

call +44 (0) 1274 851 999
email sales.uk@ferno.com
www ferno.com



FERNO Compact 2 Track In-vehicle storage bracket




Secured, ready for use and easy to deploy

If you carry the Compact 2 Track in your Ambulance this crash-tested two-piece storage bracket is an essential extra.

The track fits in first and can stay in the bracket if you need to use the Compact 2 purely as a carry chair, the chair locks in over the track and is secured further with a strap. CEN Crash tested to 10G and CEN compliant, you can rest assured that when stored, the equipment is safe and secured to the highest standard.

- **Purpose** – designed to secure the Compact 2 Track in vehicle
- **Protects** – tested and certified to BS EN 1789
- **Easy to use** – the track and chair simply click into place and are easily deployed
- **Easy to clean and maintain** – ask about our comprehensive service contract for all equipment

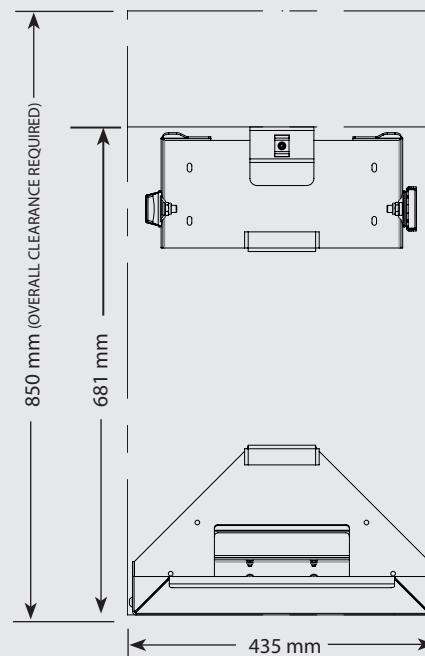

Part Nos:

009103800

Compact 2 Track 2-Piece In-Vehicle Storage Bracket

FERNO - C2T Bracket FPC010 03.25

call +44 (0) 1274 851 999
email sales.uk@ferno.com
www ferno.com



Ferno (UK) Limited, Ferno House, Stubbs Beck Lane, Cleckheaton, West Yorkshire, BD19 4TZ, England