



# UNI-ClipDoc System™



## Developed to create seamless patient transfers between Ground and Air Transport

Rescue professionals who deal with the transfer of patients will appreciate the importance of making any vehicle transition as simple and stress-free as possible.

Between air and ground, FERNO's UNI-ClipDoc is the solution.


This highly engineered interface allows equipment such as the LifePort AeroSled TS platform and Incubator interfaces to be seamlessly transferred from Aircraft to Ambulance trolley and safely transport the patient by road ambulance for the remainder of the journey.

The UNI-ClipDoc simply mirrors the Aircraft/Helicopter floor, but as you would expect from FERNO the UNI-ClipDoc has been rigorously tested and has passed all required safety standards - BS EN 1789.

Compatible with all Ferno's range of ambulance trolleys, Uni-ClipDoc gives complete interchangeability across various transportation platforms.

- **Lightweight yet incredibly strong** – manufactured from aluminium and stainless steel
- **Universal and adaptable** – a Ratchet Clamp System locks UNI-ClipDoc to traditional tubular-framed trolleys such as Pegasus, Megabus and Harrier, a Quick Release Clamp system attaches UNI-ClipDoc to the side interface rails on the latest generation of trolleys such as iNX and POWERX
- **Saves vital rescue and transfer time** – now a patient can remain safe and secure in the equipment they were originally transported in, eliminating the need to carry out transfers from one platform to another, reducing time and potential manual handling risks.

# UNI-ClipDoc System™

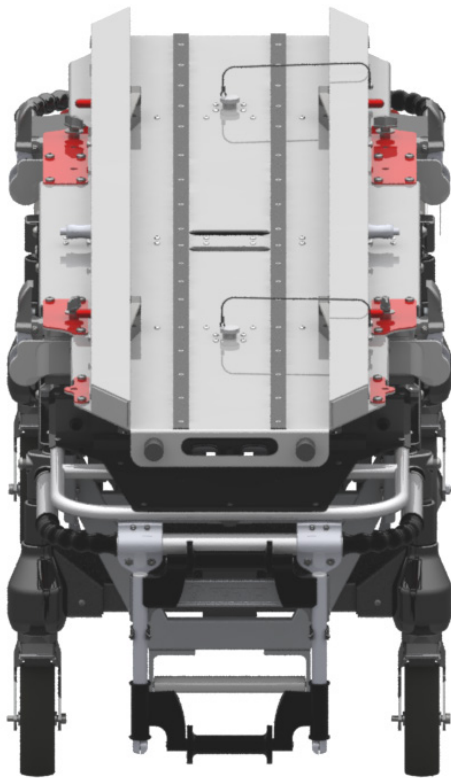


## Features & Benefits

- **Reduction in manual handling injuries** – less lifting during transfer also ensures the patient remains stable
- **Enhanced safety** – the UNI-ClipDoc is CEN compliant, dynamically tested and certified to BS EN 1789 on the FERNO Ambulance Trolley range
- **Custom-made** – Each UNI-ClipDoc will require an individual technical assessment to determine pin centres and whether a bridging piece is needed between trolley and craft to aid loading/unloading
- **Quick, easy assembly** – and for ease of storage the UNI-ClipDoc comes in two halves with sturdy carry bags.

## Specifications

<b>LENGTH</b>	1955 mm
<b>WIDTH</b>	585 mm
<b>DEPTH</b>	205 mm
<b>WEIGHT</b>	32.3 kg
<b>FULL KIT FOLDED in CARRY BAG</b>	Approx. 1200 x 800 x 210 mm, 33 kg

## Accessories

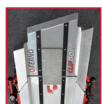

- Load/Unload extension platform made to fit your aircraft's specific requirements

## Part Numbers

031379000 UNI-ClipDoc base kit (including fixation brackets and straps for all Ferno's range of trolleys)  
TBC Load/unload extension platform

