

## TRANSCEND™

### EVAC CHAIR



The **Transcend™ Evac Chair** helps you quickly and safely transport a seated passenger down stairs and over flat surfaces in an evacuation situation.

Belted tracks enable operators to glide the chair down stairs instead of carrying it.

The durable Transcend Evacuation Chair can carry up to 550 lb (250 kg).



#### EZ GLIDE®

technology provide a comfortable ride on almost any surface with 60% less shock than competitive stair chairs.



#### Extending Footrest

provides added support and comfort for patients.




#### Rear Lift Handles

with improved ergonomics for lifting or carrying the chair

- 15-year warranty

# TRANSCEND™

## EVAC CHAIR

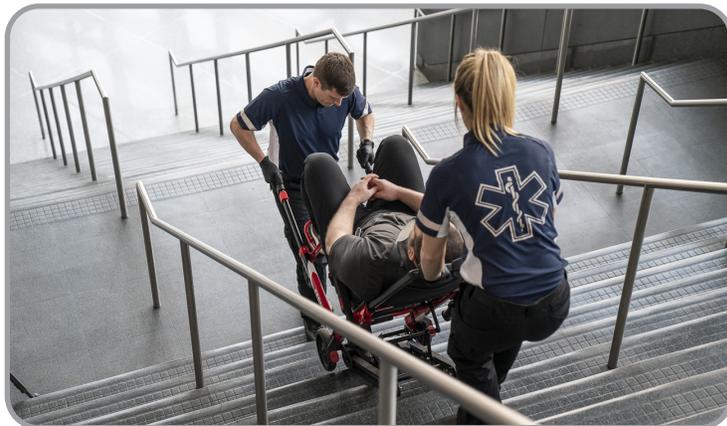

## FEATURES AND BENEFITS

- **All-metal frame construction** provides durability and long life
- **Smooth 5-inch front swivel wheels and extra tall 8-inch rear wheels** with integrated shock absorbing tires that navigate rough terrain and do not collect debris
- **One inch track-to-ground clearance** provides smooth rolling over carpet and rough surfaces, as well as easier descent down stairs
- **5-position lift bar and handles** reduce the strain of lifting and provide maximum control for operators of various heights
- **Extending footrest and lower track angle** increase patient sense of security
- **Patient seating surface** is contoured and has removable injection molded seat panels
- **Positive locking hinge mechanism** ensures chair is locked in place whether deployed for use or folded for carrying and storage
- **15-year warranty**

## SPECIFICATIONS

<b>Item Number</b>	Base model: 0731404
<b>Height (min)</b>	40 in / 1022 mm
<b>Height (max)</b>	64 in / 1633 mm
<b>Width (overall)</b>	20 in / 505 mm
<b>Width (seat)</b>	16 in / 408 mm
<b>Depth (front to back)</b>	30 in / 760 mm
<b>Length (folded)</b>	40 in / 912 mm
<b>Width (folded)</b>	20 in / 507 mm
<b>Depth (folded)</b>	7.9 in / 200 mm
<b>Weight</b>	37 lb / 16.8 kg
<b>Load Limit</b>	550 lb / 250 kg

## TRANSCEND OPTIONS



Wheel Locks



Armrests



Extending Lift Bar



Footrest

## PATIENT HANDLING ACCESSORIES



Headrest  
No. 0822023



Linkage/Close Handle  
Kit No. 0822536



Rear Lift Handles  
No. 0822028



Restraint Set  
No. 0822024